




ADVACRYL RAPID

Fast Absorbing Polyglactin-910 (Undyed)

Braided Coated Synthetic Absorbable Surgical Suture

Accelerated performance for confidence
in delicate tissue approximation



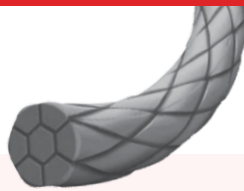

Accelerated performance for confidence in delicate tissue approximation

Product Description

Sr. No.	Code	Suture			Needle				Units Per Box
		Colour	Size	Length (cm)	Curvature	Type	Length (mm)	Speciality	
1	ADVA2721	UNDYED	0	60	3/8 CIRCLE		35	-	12
2	ADVA2763	UNDYED	0	110	1/2 CIRCLE		36	-	12
3	ADVA2763FS	UNDYED	0	100	1/2 CIRCLE		36	-	12
4	ADVA2763S	UNDYED	0	70	1/2 CIRCLE		36	-	12
5	ADVA2763SP	UNDYED	0	100	1/2 CIRCLE		40	-	12
6	ADVA2763RB	UNDYED	0	100	1/2 CIRCLE		36	LONG TAPER	12
7	ADVA9963	UNDYED	0	90	1/2 CIRCLE		36	-	12
8	ADVA9972	UNDYED	0	90	1/2 CIRCLE		40	HEAVY	12
9	ADVA9980	UNDYED	0	90	1/2 CIRCLE		48	-	12
10	ADVA9987	UNDYED	0	90	1/2 CIRCLE		40	-	12
11	ADVA9964	UNDYED	1	90	1/2 CIRCLE		36	-	12
12	ADVA2720	UNDYED	2-0	60	3/8 CIRCLE		36	-	12
13	ADVA2761	UNDYED	2-0	100	1/2 CIRCLE		36	-	12
14	ADVA2761SP	UNDYED	2-0	100	1/2 CIRCLE		40	-	12
15	ADVA2762	UNDYED	2-0	90	1/2 CIRCLE		35	-	12
16	ADVA2762S	UNDYED	2-0	30	1/2 CIRCLE		35	-	12
17	ADVA2762ER	UNDYED	2-0	110	1/2 CIRCLE		35	-	12
18	ADVA9933	UNDYED	2-0	75	3/8 CIRCLE		26	PRIME	12
19	ADVA9938	UNDYED	2-0	75	3/8 CIRCLE		36	PRIME	12
20	ADVA9962	UNDYED	2-0	90	1/2 CIRCLE		36	-	12
21	ADVA9970	UNDYED	2-0	70	3/8 CIRCLE		24	-	12
22	ADVA9973	UNDYED	2-0	90	1/2 CIRCLE		37	HEAVY	12
23	ADVA9974	UNDYED	2-0	70	1/2 CIRCLE		26	-	12
24	ADVA9975CR	UNDYED	2-0	75	1/2 CIRCLE		31	LONG TAPER	12
25	ADVA9981	UNDYED	2-0	75	1/2 CIRCLE		30	-	12
26	ADVA9982	UNDYED	2-0	75	1/2 CIRCLE		35	-	12
27	ADVA9983	UNDYED	2-0	75	1/2 CIRCLE		45	-	12
28	ADVA9988	UNDYED	2-0	90	1/2 CIRCLE		48	-	12
29	ADVA2719	UNDYED	3-0	45	3/8 CIRCLE		16	-	12
30	ADVA2732	UNDYED	3-0	70	3/8 CIRCLE		26	-	12
31	ADVA2732S	UNDYED	3-0	45	3/8 CIRCLE		26	-	12
32	ADVA2735	UNDYED	3-0	70	1/2 CIRCLE		22	-	12
33	ADVA2764	UNDYED	3-0	70	1/2 CIRCLE		26	-	12
34	ADVA9919	UNDYED	3-0	75	3/8 CIRCLE		16	PRIME	12
35	ADVA9931	UNDYED	3-0	75	3/8 CIRCLE		19	-	12
36	ADVA9932	UNDYED	3-0	75	3/8 CIRCLE		26	PRIME	12
37	ADVA9935	UNDYED	3-0	45	1/2 CIRCLE		22	-	12
38	ADVA9974CR	UNDYED	3-0	75	1/2 CIRCLE		22	LONG TAPER	12
39	ADVA9975	UNDYED	3-0	70	3/8 CIRCLE		19	-	12
40	ADVA9976	UNDYED	3-0	70	3/8 CIRCLE		24	-	12
41	ADVA9984	UNDYED	3-0	45	1/2 CIRCLE		22	PRIME	12
42	ADVA9986	UNDYED	3-0	90	3/8 CIRCLE		26	-	12

Round Bodied

Taper Cut

Reverse Cutting

Cutting

Sr. No.	Code	Suture			Needle				Units Per Box
		Colour	Size	Length (cm)	Curvature	Type	Length (mm)	Speciality	
43	ADVA9991	UNDYED	3-0	45	3/8 CIRCLE		19	-	12
44	ADVA9993	UNDYED	3-0	75	3/8 CIRCLE		24	-	12
45	ADVA9995	UNDYED	3-0	75	3/8 CIRCLE		19	PRIME	12
46	ADVA2718	UNDYED	4-0	45	3/8 CIRCLE		16	PRIME	12
47	ADVA2718RB	UNDYED	4-0	45	3/8 CIRCLE		16	LONG TAPER	12
48	ADVA2718RC	UNDYED	4-0	45	COMPOUND CURVE		16	-	12
49	ADVA9918	UNDYED	4-0	75	3/8 CIRCLE		16	PRIME	12
50	ADVA9930	UNDYED	4-0	75	3/8 CIRCLE		19	-	12
51	ADVA9977	UNDYED	4-0	70	3/8 CIRCLE		16	-	12
52	ADVA9978	UNDYED	4-0	70	1/2 CIRCLE		17	-	12
53	ADVA9985	UNDYED	4-0	75	1/2 CIRCLE		19	PRIME	12
54	ADVA9989	UNDYED	4-0	45	1/2 CIRCLE		19	PRIME	12
55	ADVA9915	UNDYED	5-0	45	3/8 CIRCLE		11	PRIME	12
56	ADVA9916	UNDYED	5-0	45	1/2 CIRCLE		12	PRIME	12
57	ADVA9969	UNDYED	5-0	75	1/2 CIRCLE		17	LONG TAPER	12
58	ADVA9979	UNDYED	5-0	70	3/8 CIRCLE		12	-	12
59	ADVA9979S	UNDYED	5-0	45	3/8 CIRCLE		21	-	12
60	ADVA9990	UNDYED	5-0	45	1/2 CIRCLE		16	LONG TAPER	12
61	ADVA9992	UNDYED	5-0	75	3/8 CIRCLE		12	-	12
62	ADVA9994	UNDYED	6-0	45	3/8 CIRCLE		11	PRIME	12

Sr. No.	Code	Suture			Needle				Units Per Box
		Colour	Size	Length (cm)	Curvature	Type 1	Type 2	Length (mm)	
63	ADVA2763DA	UNDYED	0	70	1/2 CIRCLE			36	12
64	ADVA9963DA	UNDYED	0	36	1/2 CIRCLE			36	12
65	ADVA9971	UNDYED	0	90	1/2 CIRCLE			43	12
66	ADVA2666	UNDYED	2-0	90	1/2 CIRCLE			40	12
67	ADVA2777	UNDYED	2-0	140	1/2 CIRCLE			36	12
68	ADVA2777TR	UNDYED	2-0	140	1/2 CIRCLE			36 & 40	12
69	ADVA9979S	UNDYED	5-0	45	3/8 CIRCLE			21	12

Round Bodied

Taper Cut

Reverse Cutting

Cutting

Product Description


Composition	90% Glycolide, 10% L-lactide
Coating	Copolymer of Poly (glycolide-co-lactide) and Calcium Stearate
Size	Size 1 to 6-0
Colour	Undyed
Absorption Profile	Essentially complete up to 42 days
Needle Availability	A wide variety of Round Bodied (Taper Point), Ball Point, Cutting, Reverse Cutting, Taper Cut, Trocar, Blunt Point, Oval Round Bodied, Diamond Point & Spatula Micro Point needles

Tensile Strength¹

Post implantation approximately 42% of the original tensile strength remains at 5 days. All of the original tensile strength is lost by approximately 10 to 14 days post implantation.

The characteristic rapid loss of strength is achieved by use of a polymer material with a lower molecular weight than regular ADVACRYL (polyglactin 910) suture.

Knot pull tensile strength (KgF)²

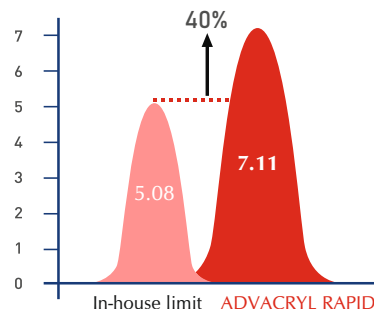

REFERENCE _____


1 : Data on file with AMS

2 : Data on file with AMS - DV done on Polyglactin 910, size - 1, 36 mm Taper cut needle

Multifilament Material




Undyed

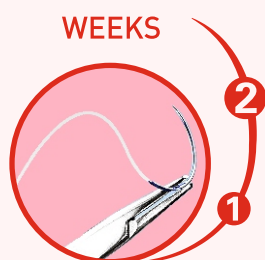

Provides excellent support for a range of critical wound healing period to secure closure in a variety of tissues, procedures and patients.

Premium Needles for consistent Superior Performance

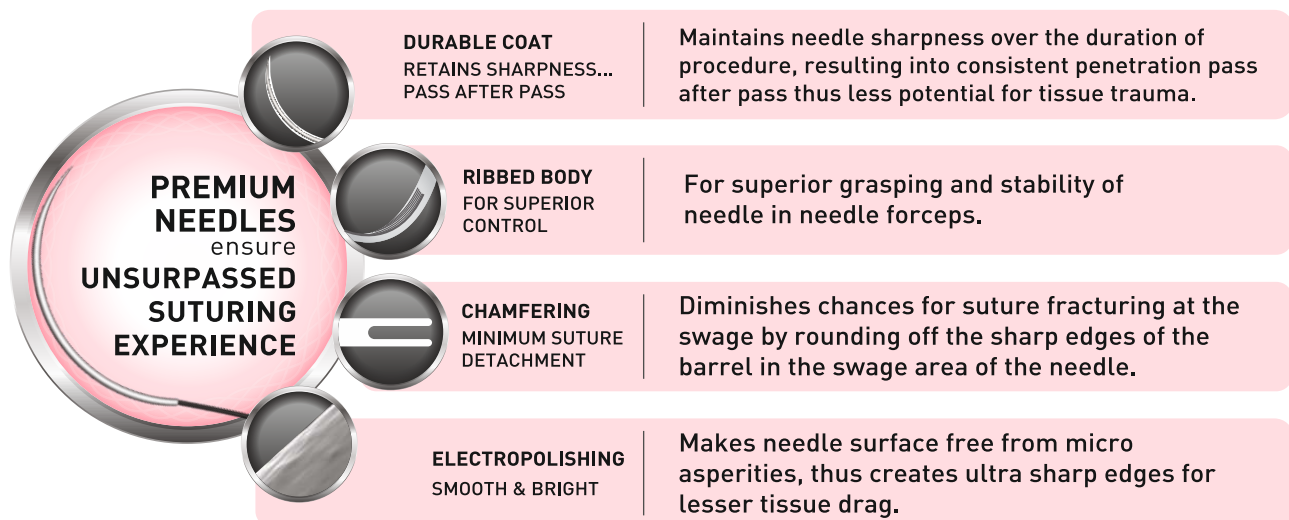

Minimal tissue trauma enables faster recovery



P-21-22, 25-26 & 34-35, GIDC, Manjusar, Tal - Savli, Dist. Vadodara - 391775, Gujarat, India
T : (91) 2667 671302 | amsltd.com | info@amsltd.com

Singapore Office : 111, Somerset Road, 07-08, 111 Somerset, Singapore - 238164.
T: (65) 636910 80 / 89 | E: info@amsltd.com | amsltd.com